



CZR PRESS TRIAL-L KIT

CZR PRESS TRIAL-L KIT is a necessary kit for using INGOT-L, including SHADE BASE STAIN, INGOT-L, CZR PORCELAIN, and PLUNGER. It contains 5 shade colors, which are A1, A2, A3, NW0, and NW0.5. The beginning user of CZR PRESS should also prepare for NORITAKE PRESS INVESTMENT, RING/FORMER, IS LIQUID, and CZR FORMING LIQUID.

L-Ingot
for use when utilizing the "Layering Method" & "LF Layering Method"

INGOT-L	10g each	5 colors
SB STAIN	3g each	5 colors
CZR PORCELAIN POWDER	10g each	8 colors
PLUNGER		1 set
SB STAIN COLOR GUIDE		1 piece
CZR PORCELAIN INSTRUCTION		1 copy
TECHNICAL INSTRUCTIONS		1 copy



CZR PRESS TRIAL-H KIT

CZR PRESS TRIAL-H KIT is a necessary kit for using INGOT-H, including SHADE BASE STAIN, INGOT-H, CZR PRESS GLAZE, CZR ES (A+), and PLUNGER. It contains 5 shade colors, which are A1, A2, A3, NW0, and NW0.5. The beginning user of CZR PRESS should also prepare for NORITAKE PRESS INVESTMENT, RING/FORMER, IS LIQUID, and ES LIQUID.

H-Ingot
for use when utilizing the "Staining Method" & "LF Layering Method"

INGOT-H	10g each	5 colors
SB STAIN	3g each	5 colors
CZR ES A+	3g	1 color
CZR PRESS GLAZE	10g	1 color
PLUNGER		1 set
SB STAIN COLOR GUIDE		1 piece
CZR PORCELAIN INSTRUCTION		1 copy
TECHNICAL INSTRUCTIONS		1 copy



CZR PRESS INGOT-L KIT

CZR PRESS INGOT-L KIT is a low translucent type of press ceramic for zirconia framework. It is developed to recreate the margin area precisely by pressing INGOT-L. INGOT-L is a press ceramic using with the layering technique. After INGOT-L is pressed, the crown can be made by using CZR TRANSLUCENT PORCELAIN or CZR LUSTER PORCELAIN. Can make a single anterior crown without frameworks using with CZR PRESS LF porcelain.

BASIC COLOR	10g each	20 colors
PLUNGER		1 set
TECHNICAL INSTRUCTIONS		1 copy

5g Ingot is available as individual item.



CZR PRESS INGOT-H KIT

CZR PRESS INGOT-H KIT is a high translucent type of press ceramic for zirconia framework. It is developed to recreate the margin area precisely by pressing INGOT-H. INGOT-H is a press ceramic using with the stain technique. After INGOT-H is pressed, the crown can be made by using CZR EXTERNAL STAIN or CZR PRESS GLAZE POWDER. CZR PRESS LF porcelain can be used after pressing also.

BASIC COLOR	10g each	20 colors
PLUNGER		1 set
TECHNICAL INSTRUCTIONS		1 copy

5g Ingot is available as individual item.



CZR PRESS SHADE BASE STAIN KIT

CZR PRESS SHADE BASE STAIN KIT is a kind of the stain to make the basic shade color of zirconia framework before pressing CZR PRESS INGOT. The basic shade color can be recreated thinly, transparently and naturally with the simple stain technique.

BASIC COLOR	3g each	20 colors
IS LIQUID	10ml	1 bottle
SB STAIN COLOR GUIDE		1 piece
TECHNICAL INSTRUCTIONS		1 copy



CZR PRESS ESTHETIC WHITE KIT

CZR PRESS ESTHETIC WHITE KIT is developed to create whiter shades rather than NW shades with Press technique. It is not only pure bleaching shades but also natural and beautiful white shades.

SHADE BASE STAIN	3g each	SS White, SSB1
ESTHETIC WHITE INGOT	10g each	EW00, EW0, EW, EWY
EXTERNAL STAIN	3g	B+
GLAZE	10g	CZR Press Glaze
TECHNICAL INSTRUCTION		1 copy



CZR PRESS LAYERING KIT


CZR PRESS LAYERING KIT is a set of CZR porcelains selected for using CZR PRESS INGOT with zirconia frameworks.

BODY	10g each	20 colors
ENAMEL		3 colors
TRANSLUSENT		3 colors
MODIFIER		14 colors
ADD-ON		2 colors
LUSTER		8 colors
MRP		1 color
COLOR TABLE/WITHOUT SB		1 piece
TECHNICAL INSTRUCTIONS		1 copy



CZR PRESS GLAZE POWDER

This CZR PRESS GLAZE POWDER is necessary to produce glossiness of the crown in the H-ingots staining technique. Be sure to use this glaze powder and refer to our CZR PRESS technical instruction for the glazing temperature. This glaze powder is available individually. **The use of IS Liquid is highly recommended.**

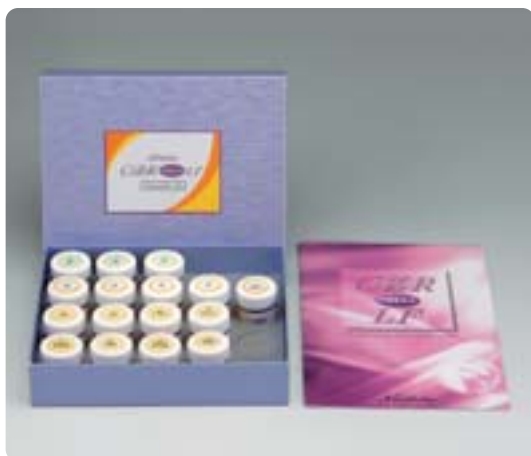

CZR PRESS LF ASSORTMENT KIT

Noritake CZR PRESS LF is low fusing porcelain to be built up to make an enamel layer after pressing CZR PRESS ingot. By using this porcelain with pressed ceramic without a Zirconia framework, you can make a single crown, a porcelain laminate veneer, an inlay and an onlay. Please use Noritake LF Liquid for this porcelain.

ENAMEL	10g each	3 colors
TRANSLUSENT/LUSTER		12 colors
AD-T		1 color
TECHNICAL INSTRUCTIONS		1 copy

CZR PRESS LF INTERNAL STAIN KIT



CZR PRESS LF INTERNAL STAIN KIT is a set of Low fusing stains for CZR PRESS LF porcelain. It is necessary to use this stain before CZR PRESS LF layering. Baking temperature of this stain is 800°C.

1. WHITE	7. CERVICAL 2	13. A+	3g each
2. INCISAL BLUE 1	8. CERVICAL 3	14. B+	
3. INCISAL BLUE 2	9. EARTH BROWN	15. C+	
4. MAMELON ORANGE 1	10. REDDISH BROWN	16. D+	
5. MAMELON ORANGE 2	11. SALMON PINK		
6. CERVICAL 1	12. RED		
BRIGHT (DILUTION)	3g	1 color	
IS LIQUID	10ml	1 bottle	
IS COLOR GUIDE		1 piece	
TECHNICAL INSTRUCTIONS		1 copy	

CZR PRESS LF EXTERNAL STAIN KIT



CZR PRESS LF EXTERNAL STAIN KIT is a set of Low fusing stains for CZR PRESS LF porcelain. It is necessary to use this stain after CZR PRESS LF layering. Baking temperature of this stain is 840°C.

1. PURE WHITE	8. ORANGE 1	15. SALMON PINK	3g each
2. GRAY	9. ORANGE 2	16. PINK	
3. BLACK	10. CERVICAL 1	17. RED	
4. BLUE	11. CERVICAL 2	18. A+	
5. GREEN 1	12. CERVICAL 3	19. B+	
6. GREEN 2	13. EARTH BROWN	20. C+	
7. YELLOW	14. REDDISH BROWN	21. D+	
GLAZE	10g	1 color	
ES LIQUID	10ml	1 bottle	
ES COLOR GUIDE		1 piece	
TECHNICAL INSTRUCTIONS		1 copy	

LF LIQUID

CZR PRESS LF porcelain should EXCLUSIVELY be mixed with LF LIQUID.

100ml

RING / FORMER SET



NORITAKE RING/FORMER SET is used for investing waxed-up crowns. 3 sizes are available, 100g, 200g and 300g.

RING / FORMER 100 - for 2g Ingot
200 - for 2g Ingot
300 - for 5g Ingot

RING/FORMER SET 100/200

RING 100	FORMER 100
RING 200	FORMER 200

RING/FORMER SET 100

RING 100	FORMER 100
----------	------------

RING/FORMER SET 200

RING 200	FORMER 200
----------	------------

RING/FORMER SET 300

RING 300	FORMER 300
----------	------------

NORITAKE PLUNGER



NORITAKE PLUNGER is a plunger used when pressing CZR PRESS INGOT.

3 pieces

NORITAKE DISPO PLUNGER



Noritake DISPO PLUNGER is disposable plunger set. 2 kinds of size are available for each 2g and 5g CZR PRESS INGOT.

for 2g Ingot	∅12mm/H37mm	50 pcs.
for 5g Ingot	∅16.5mm/H42mm	50 pcs.